

京天成生物技术(北京)有限公司 AbMax Biotechnology Co., LTD 北京经济技术开发区科创十四街 99 号 18 号楼 2 单元 201 99 Kechuang 14th St, Bldg 18, 2-201 BDA, Beijing, China 101111 Tel: 86-10-59755729, Fax: 86-10-59755718 E-mail:info@antibodychina.com www.antibodychina.com

Anti-2019-nCoV Nucleoprotein (N), Near CT


(rabbit polyclonal lgGs) Catalog # 06-0051

Immunogen: synthetic peptide corresponding to a sequence near the middle of 2019-nCoV N protein. Specificity: Recognizes the N proteins from both 2019-nCoVand SARS. Formulation: 0.05mg of rabbit polyclonal IgG of 0.01M PB, pH 7.4 Protein was determined by OD₂₈₀ absorbance. Physical State: Frozen liquid Purity: > 95% based on SDS-PAGE **Storage:** Store the product at -20 $^{\circ}$ C. Product is stable for about 6 weeks at 2-8 $^{\circ}$ C as an undiluted liquid. Prepare working dilution fresh each day. Avoid repeated freezing and thawing.

FOR IN VITRO RESEARCH USE ONLY NOT FOR USE IN HUMANS

Research Applications

Indirect ELISA: to be used at 1~2µg/ml. **Western Blot:** to be used at 0.5~2µg/ml

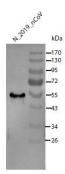

Immunohistochemistry: to be used at 10-30µg/ml.